

Enclosed Plate and Frame Filter BECO INTEGRA® PLATE A600 EC

Stainless Steel

BECO INTEGRA PLATE A600 EC enclosed plate and frame filter is the depth filtration system for critical filtration tasks in the pharmaceutical and biotech, chemical, cosmetics, and food industries.

BECO INTEGRA PLATE A600 EC is a multi-sheet-filter with an external chamber. The product channels are sealed with depth filter sheets only. The sealing of the whole filter package is realized by means of gaskets around the external chamber.

Special Features

- **CIP** of filter package without depth filter sheets.
- **Complete drainage** of filter due to unique filter plate design.
- **Safe operation:** No drip loss due to sealed external chamber and automatic hydraulic pressure regulation.
- **Easy handling** and **good cleanability** thanks to ergonomic design of chassis, and open support grid of filter plates; cake discharge downwards into tray underneath of filter package.
- **High flexibility:** The filter can be used as a sheet filter or as a pre-coat filter. For pre-coating/removing of high particle load, cake frames in different widths are available.
- BECO INTEGRA PLATE A600 EC filter is designed and manufactured according to the **cGMP guidelines** and the **FDA regulations**.

Customization or additional components allow for system expansion or complete systems according to customer specifications.



Technical Data

Operating pressure:	max. 87 psi (600 kPa/6 bar)
Differential pressure:	max. 58 psi (400 kPa/4 bar)
Operating temperature:	max. 284 °F (140 °C)
Filter area per sheet:	3.6 ft ² (0.33 m ²)

Materials

Product contacting parts:	Stainless steel AISI 316L/AISI 316Ti
Other parts:	Stainless steel AISI 304
Gaskets/O-rings:	Silicone, EPDM, Viton



Powering Business Worldwide

Type	Filter area		Cake volume*		Chassis size	Filter plates/ cake frames	Depth filter sheets 23.0 x 26.8 in (584 x 680.5 mm)	Weight	
	ft ²	(m ²)	gal	(l)				qty	qty
A600 SF	215.3	(20.00)	-		1	29/30	60	3,477	(1,577)
	387.5	(36.00)	-		2	53/54	108	4,306	(1,953)
	538.2	(50.00)	-		3	74/75	150	4,996	(2,266)
A600 ASF with 0.31" (8 mm) cake frames	215.3	(20.00)	7.9	(30)	1	29/30	60	3,448	(1,564)
	387.5	(36.00)	14.3	(54)	2	53/54	108	4,259	(1,932)
	538.2	(50.00)	19.8	(75)	3	74/75	150	4,930	(2,236)
A600 ASF with 0.71" (18 mm) cake frames	150.7	(14.00)	24.0	(91)	1	20/21	42	3,404	(1,544)
	272.7	(25.33)	43.0	(164)	2	37/38	76	4,182	(1,897)
	380.3	(35.33)	60.5	(229)	3	52/53	106	4,833	(2,192)
A600 ASF with 0.98" (25 mm) cake frames	129.2	(12.00)	31.7	(120)	1	17/18	36	3,455	(1,567)
	229.6	(21.33)	56.3	(213)	2	31/32	64	4,246	(1,926)
	315.7	(29.33)	77.4	(293)	3	43/44	88	4,892	(2,219)
A600 ASF with 1.57" (40 mm) cake frames	93.3	(8.67)	39.9	(151)	1	12/13	26	3,475	(1,576)
	165.0	(15.33)	70.8	(268)	2	22/23	46	4,268	(1,936)
	229.6	(21.33)	98.5	(373)	3	31/32	64	4,956	(2,248)

* usable cake volume

Filter Plates/Cake Frames

	Width	Cake volume*		Weight	
	in (mm)	gal	(l)	lb	(kg)
Filter plate	0.31 (8)			14.3	(6.5)
Feeding plate	0.31 (8)			14.3	(6.5)
Cake frame	0.31 (8)	0.26	(1.00)	13.5	(6.1)
Cake frame	0.71 (18)	1.14	(4.33)	23.3	(10.5)
Cake frame	0.98 (25)	1.76	(6.66)	32.2	(14.6)
Cake frame	1.57 (40)	3.08	(11.66)	51.6	(23.4)


* usable cake volume

Connections

Filter inlet:	DN 65, DIN 11851/ISO 2 ½"
Filter outlet:	DN 65, DIN 11851/ISO 2 ½"
External chamber:	DN 15, DIN 11851/ISO ¾"

Other Designs

Flange according to DIN 2633/ISO or tri-clamp on request

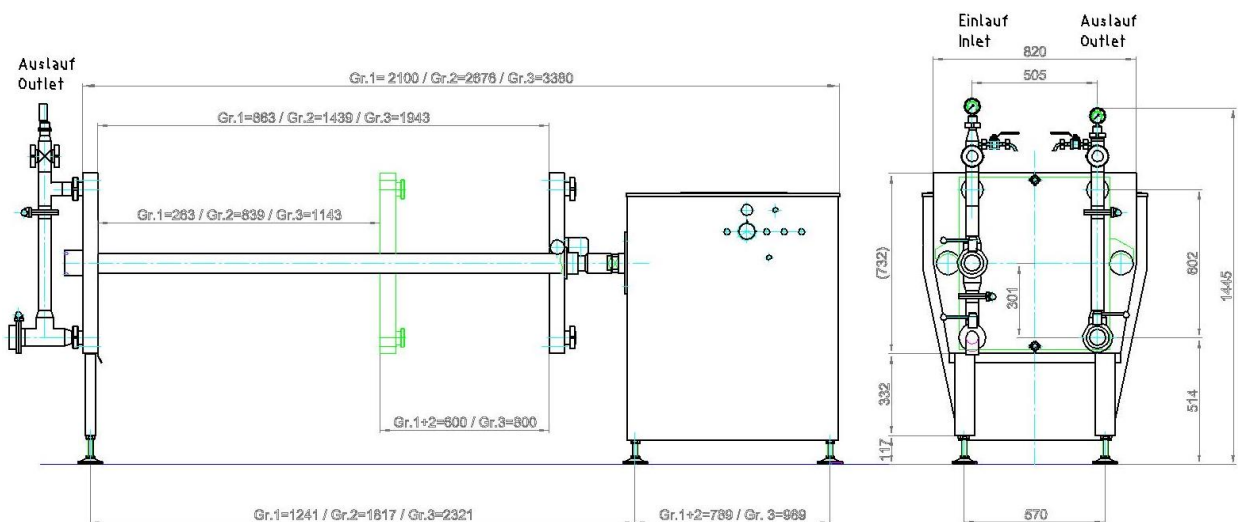

North America

44 Apple Street
Tinton Falls, NJ 07724
Toll Free: 800 656-3344
(North America only)
Tel: +1 732 212-4700

Europe/Africa/Middle East

Auf der Heide 2
53947 Nettersheim, Germany
Tel: +49 2486 809-0

Friedensstraße 41
68804 Altlußheim, Germany
Tel: +49 6205 2094-0

An den Nahewiesen 24
55450 Langenlonsheim, Germany
Tel: +49 6704 204-0

China

No. 3, Lane 280,
Linhong Road
Changning District, 200335
Shanghai, P.R. China
Tel: +86 21 5200-0099

Singapore

100G Pasir Panjang Road #07-08
Singapore 118523
Tel: +65 6825-1668

Brazil

Rua Clark, 2061 - Macuco
13279-400 - Valinhos, Brazil
Tel: +55 11 3616-8400

**For more information, please
email us at filtration@eaton.com
or visit www.eaton.com/filtration**

EN
1 A 2.6.11.2
03-2019

© 2019 Eaton. All rights reserved. All trademarks and registered trademarks are the property of their respective owners. All information and recommendations appearing in this brochure concerning the use of products described herein are based on tests believed to be reliable. However, it is the user's responsibility to determine the suitability for his own use of such products. Since the actual use by others is beyond our control, no guarantee, expressed or implied, is made by Eaton as to the effects of such use or the results to be obtained. Eaton assumes no liability arising out of the use by others of such products. Nor is the information herein to be construed as absolutely complete, since additional information may be necessary or desirable when particular or exceptional conditions or circumstances exist or because of applicable laws or government regulations.



Powering Business Worldwide